


Background

Kras is a member of the RAS protein family, which are a class of small GTPases involved in cell signaling pathways. The Ras signaling pathway regulates diverse cellular processes, including cell proliferation, differentiation, and survival. Conversion of Ras from the inactive GDP-bound state to the active GTP-bound state activates the downstream effector and promotes cell growth. RAF is a key downstream effector of RAS. Since the frequently mutated *Ras* genes are associated with various human tumors, the Ras-RAF signaling pathway is considered an important therapeutic target for cancer treatment. However, Ras is considered undruggable since it lacks suitable binding pockets on the surface. Recently, a discovery of a small molecule inhibitor blocks Ras-RAF signaling pathway by remodeling Cyclophilin A (CYPA) and forming a CYPA:drug:KRAS ternary complex. This inhibitory strategy provides a new method for developing drugs targeting Kras for treatment of cancers.

Assay Principle

The Kras (G12C) Inhibitor assay kit is a TR-FRET based assay, which is designed to screen Kras inhibitors and determine the Kras-inhibitor binding affinity. Tag2-Kras (G12C) in this assay kit is loaded with GppNHp, which represents the activated Kras. The Ras binding domain (RBD) of cRAF has a Tag1 at N-terminus. A Terbium-labeled anti-Tag2 antibody binding to the Tag2-Kras serves as a fluorescence donor (HTRF donor), activation of which results in fluorescence resonance energy transfer (FRET) if Tag1-cRAF binds to the Kras, since the binding brings Terbium on the anti-Tag2 antibody close to the fluorophore on the anti-Tag1 antibody (HTRF acceptor). Thus, the binding status can be quantitatively measured by calculating the ratio of the emission fluorescence intensity of the acceptor (665 nm) and donor (620 nm). If an inhibitor associated with CYPA binds to the Kras and blocks the cRAF binding, the HTRF signal will be reduced.

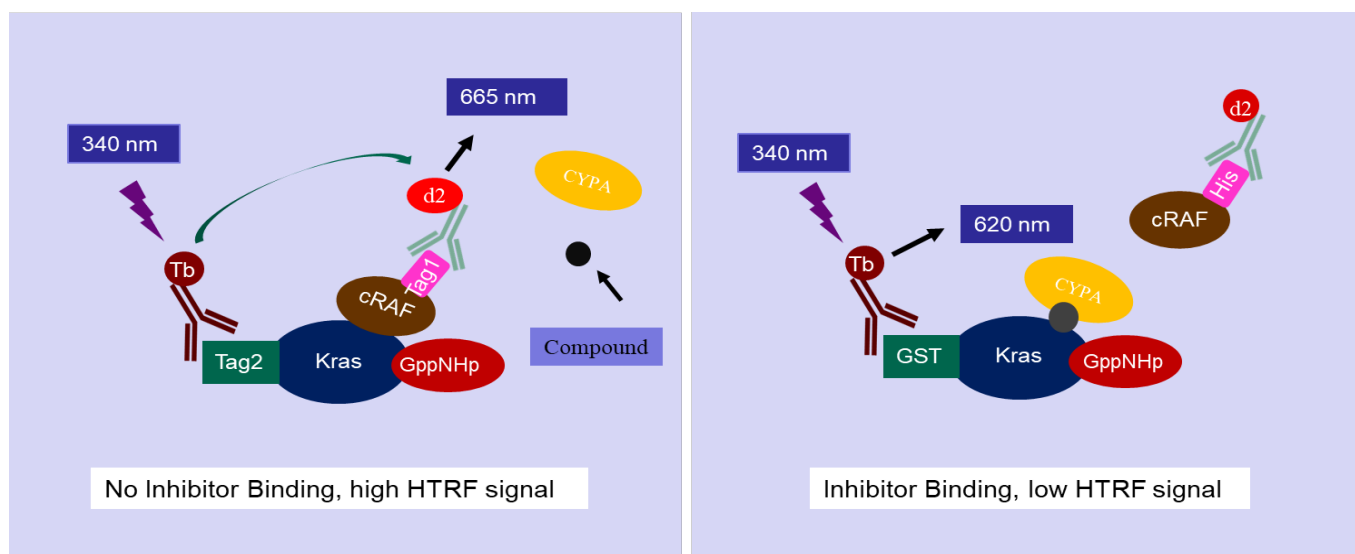

Application

High throughput screening of compounds that inhibit the binding between activated Kras (G12C) and cRAF for drug discovery.

Plate Reader

A HTRF® certified microplate reader capable of measuring Time Resolved Fluorescence Resonance Energy Transfer (TR-FRET) is required.

Components

Catalog number	Item	Amount	Storage
5727-CK-B	2X Kras Binding buffer	25 mL	-20°C
5727-4122-T2P	Recombinant human Tag2-Kras (G12C), GppNHp loaded	15 µL	-80°C
5727-2972	Recombinant Cyclophilin A (CYPA)	200 µL	-80°C
7237231-T1	Recombinant human Tag1-cRAF, RBD	8 µL	-80°C
37882	Terbium-labeled anti-Tag2 antibody	20 µL	-80°C
44732	Fluorescence labeled anti-Tag1 antibody	20 µL	-80°C
	384-well microplate	1	Room temperature

Materials needed but not supplied

1. Microplate reader, HTRF® certified microplate reader
2. 0.5 M DTT
3. Adjustable micro-pipettor
4. Sterile Tips

Assay protocol

1. Prepare compound dilution buffer containing 2 mM DTT (CD buffer)

For example, mix 250 µl of 2X Kras Binding Buffer, 179 µl of distilled water, 70 µl of CYPA and 1 µl of 0.5 M DTT. Make only enough CD buffer as needed for the assay. Store the remaining Binding buffer at -20°C.

2. Prepare the inhibitor compound solution

If the inhibitor compound is dissolved in water, make a solution of the compound 10-fold higher than the final concentration in CD buffer (since you will add 2 µl to the 20 µl reaction).

If the inhibitor compound is dissolved in DMSO, make a 100-fold higher concentration of the compound than the highest concentration you want to test in DMSO. Then make a 10-fold dilution in CD buffer (at this step, the compound concentration is 10-fold higher than the final concentration and the DMSO concentration is 10%). To determine an IC₅₀ or to test lower concentrations of the compound, prepare a series of further dilutions in CD buffer containing 10% DMSO (the final concentration of the DMSO will be 1% in all samples).

3. Prepare 1X Assay Buffer containing 2 mM DTT (AB buffer)

For example, mix 500 µl of 2X Kras Binding Buffer, 498 µl of distilled water and 2 µl of 0.5 M DTT. Make only enough AB buffer as needed for the assay. Store the remaining Binding buffer at -20°C.

4. Prepare Kras (G12C) solution

Thaw Kras protein on ice. Upon first thaw, briefly spin tube to recover the full contents at the bottom of the tube. Make aliquots of the enzyme for single use. Store remaining undiluted enzyme at -80°C.

Note: Kras protein is sensitive to freeze/thaw cycles. Limit number freeze-thaw cycles for best results. Do not re-use the diluted protein.

Dilute the Kras protein 110-fold (1 µL Kras G12C + 109 µL AB buffer).

Add 4 µl of diluted protein solution to each positive control well and inhibitor test well.

Add 4 µl of AB buffer to each of negative control well.

5. Add inhibitor

Add 2 µl of diluted compound solution to each inhibitor test well.

Add 2 µl of CD buffer to each of negative and positive control well.

6. Prepare cRAF solution

Thaw cRAF protein on ice. Upon first thaw, briefly spin tube to recover the full contents at the bottom of the tube. Make aliquots of the enzyme for single use. Store remaining undiluted protein at -80°C.

Note: cRAF protein is sensitive to freeze/thaw cycles. Limit number freeze-thaw cycles for best results. Do not re-use the diluted protein.

Dilute the cRAF protein 400-fold (1 µL cRAF + 399 µL of **AB buffer**).

Add 4 µl of diluted protein solution to each well.

7. Prepare dye solution

Dilute Terbium-labeled anti-Tag2 antibody and fluorescence-labeled anti-Tag1 antibody 1:200 in **AB buffer**. For example: 1 µl of Terbium-labeled anti-Tag2 antibody + 1 µl of fluorescence-labeled anti-Tag1 antibody + 198 µl of **AB buffer**.

Add 10 µl of this dye mixture to each well.

8. Incubate the reaction at room temperature for 30 minutes.

9. Measure fluorescent intensity

HTRF compatible microplate reader is needed to measure fluorescent intensity of the samples.

Fluorescent intensity should be measured twice:

1. Excitation wavelength at 340 nm and emission at 620 nm.
2. Excitation wavelength at 340 nm and emission at 665 nm.

Protocol Summary

Component	Negative Control	Positive Control	Inhibitor Test
AB buffer	4 µl		
Kras (G12C) protein		4 µl	4 µl
CD Buffer	2 µl	2 µl	
Inhibitor solution			2 µl
cRAF protein	4 µl	4 µl	4 µl
Dye solution	10 µl	10 µl	10 µl
Total Volume	20 µl	20 µl	20 µl

Incubate at room temperature for 30 minutes.

Data Analysis

1. Calculate sample HTRF signal of each well.


$$HTRF = \frac{\text{Fluorescent intensity at 665 nm}}{\text{Fluorescent intensity at 620 nm}} \times 10,000$$

2. Calculate percentage activity

In the absence of the compound (positive control), the sample signal (P) is defined as 100% activity. In the absence of enzyme (negative control), the sample signal (N) is defined as 0% activity. The percent activity in the presence of each compound is calculated according to the following equation: % activity = (S-N)/(P-N) X100, where S= the sample signal in the presence of the compound.

$$\% \text{ Activity} = \frac{S - N}{P - N} \times 100$$

Assay result

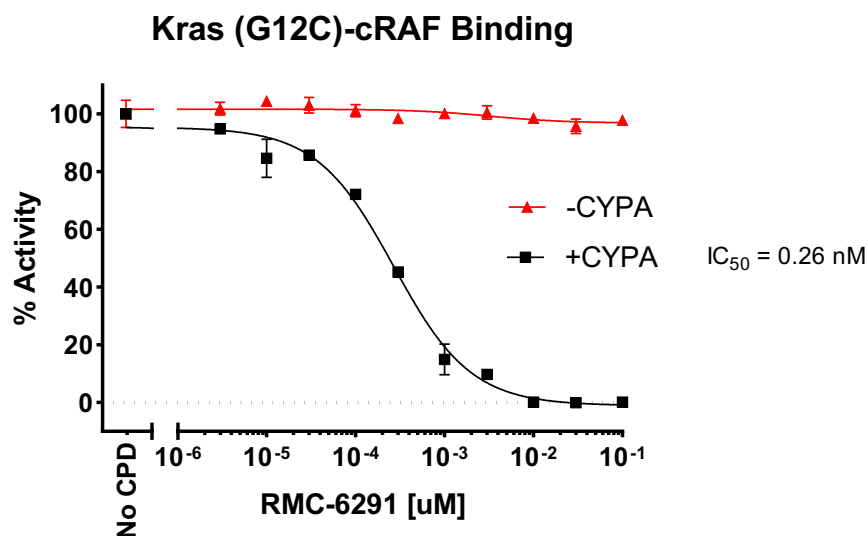

Related products:

<u>Product Name</u>	<u>Catalog #</u>	<u>Size</u>
Kras Wild Type (WT), GST-tag	5727-4121G	50 µg, 100 µg
Kras WT, GST-tag, GDP Loaded	5727-WTG-G	50 µg, 100 µg
Kras WT, GST-tag, GppNHp loaded	5727-WTG-GP	50 µg, 100 µg
Kras G12C, His -tag	5727-4122H	50 µg, 100 µg
Kras G12C, GST-tag	5727-4122G	50 µg, 100 µg
Kras G12C, GST-tag, GDP Loaded	5727-4122G -G	50 µg, 100 µg
Kras G12C, GST-tag, GppNHp loaded	5727-4122G -GP	50 µg, 100 µg
Kras G12D, GST-tag	5727-4123G	50 µg, 100 µg

Kras G12D, GST-tag, GDP Loaded	5727-4123G -G	50 µg, 100 µg
Kras G12D, GST-tag, GppNHp loaded	5727-4123G -GP	50 µg, 100 µg
Kras G12R, GST-tag,	5727-4127G	50 µg, 100 µg
Kras G12R, GST-tag, GDP Loaded	5727-4127G -G	50 µg, 100 µg
Kras G12R, GST-tag, GppNHp loaded	5727-4127G -GP	50 µg, 100 µg
Kras G12V, GST-tag,	5727-4128G	50 µg, 100 µg
Kras G12V, GST-tag, GDP Loaded	5727-4128G -G	50 µg, 100 µg
Kras G12V, GST-tag, GppNHp loaded	5727-4128G -GP	50 µg, 100 µg
Human RBD-RAF1, N-His tag, C-FLAG tag	7237231	50 µg, 100 µg
Human SOS1, No tag	7671	50 µg, 100 µg
Human SOS1, Avi-His tag	7671HA	50 µg, 100 µg
Kras WT Nucleotide Exchange Assay Kit	5727-4121NK	384 reactions
Kras G12C Nucleotide Exchange Assay Kit	5727-4122NK	384 reactions
Kras G12C Nucleotide Exchange Assay Kit	5727-4123NK	384 reactions
Kras G12R Nucleotide Exchange Assay Kit	5727-4127NK	384 reactions
Kras G12V Nucleotide Exchange Assay Kit	5727-4128NK	384 reactions
Kras WT – cRAF Binding Assay Kit	5727-4121BK	384 reactions
Kras G12C – cRAF Binding Assay Kit	5727-4122BK	384 reactions
Kras G12C– cRAF Binding Assay Kit	5727-4123BK	384 reactions
Kras G12R– cRAF Binding Assay Kit	5727-4127BK	384 reactions
Kras G12V – cRAF Binding Assay Kit	5727-4128BK	384 reactions
Kras G12C/cRAF/CYPA/Inhibitor Assay Kit	5727-4122CK	384 reactions

Service: Co-crystal Structure Determination of All KRAS Mutants

Structure Based Design, Inc. (SBD) is a sister company of **Aurora Biolabs, LLC**. Together, we have formed a strategic partnership to support early-stage drug discovery projects.

SBD has solved numerous co-crystal structures of wild-type KRAS and its various mutants. For example, we successfully determined over 30 high-resolution co-crystal structures for **Mirati Therapeutics**, with results published in peer-reviewed literature.

With established crystallization conditions and ready-to-use protein samples, we offer **rapid co-crystal structure determination services** for KRAS proteins.

Special Offer:

We now provide **co-crystal structure determination services** for the following complexes:

- **KRAS G12C/CyPA**
- **KRAS G12D/CyPA**

For more information or to get started, please contact us at:

 info@aurorabiolabs.com or info@strbd.com